

photobiologyInOut Version 0.3.1

User Guide

Pedro J. Aphalo

July 30, 2015

Contents

1	Introduction	1
2	Examples	1

1 Introduction

```
library(photobiology)
library(photobiologygg)
library(photobiologyInOut)
library(lubridate)
```

This package defines functions for importing spectral data from different instruments (Table 1) and simulation models (Table 2).

R function	Instrument	Program
read_ooss_file()	Ocean Optics spectrometers	SpectraSuite
read_oojaz_file()	Ocean Optics Jaz	<i>instrument</i>
read_macam_file()	Macam	<i>instrument</i>
read_licor_file()	LI-COR LI-1800	PC1800 (MS-DOS)

Table 1: Functions for importing measured spectral emission data

R function	Simulation model	version
read_tuv_file()	TUV (Sasha Madronich)	version 5.0 (modified)
read_tuv_file()	TUV (Sasha Madronich)	edited output

Table 2: Functions for importing simulated spectral data from models

2 Examples

Reading a file generated by Ocean Optics' Jaz spectrometer.

```
jaz.spct <- read_oojaz_file(file = "spectrum.JazIrrad")

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
```

Photon and energy irradiances.

```
q_irrad(jaz.spct, PAR()) * 1e6 # umol m-2 s-1

##          PAR
## 47.43464
## attr(,"time.unit")
## [1] "second"
## attr(,"radiation.unit")
## [1] "photon irradiance total"


e_irrad(jaz.spct, PAR()) # W m-2

##          PAR
## 10.12685
## attr(,"time.unit")
## [1] "second"
## attr(,"radiation.unit")
## [1] "energy irradiance total"
```

Plotting the spectrum.

```
plot(jaz.spct)


## Warning: Removed 1 rows containing missing values (geom_path).
```

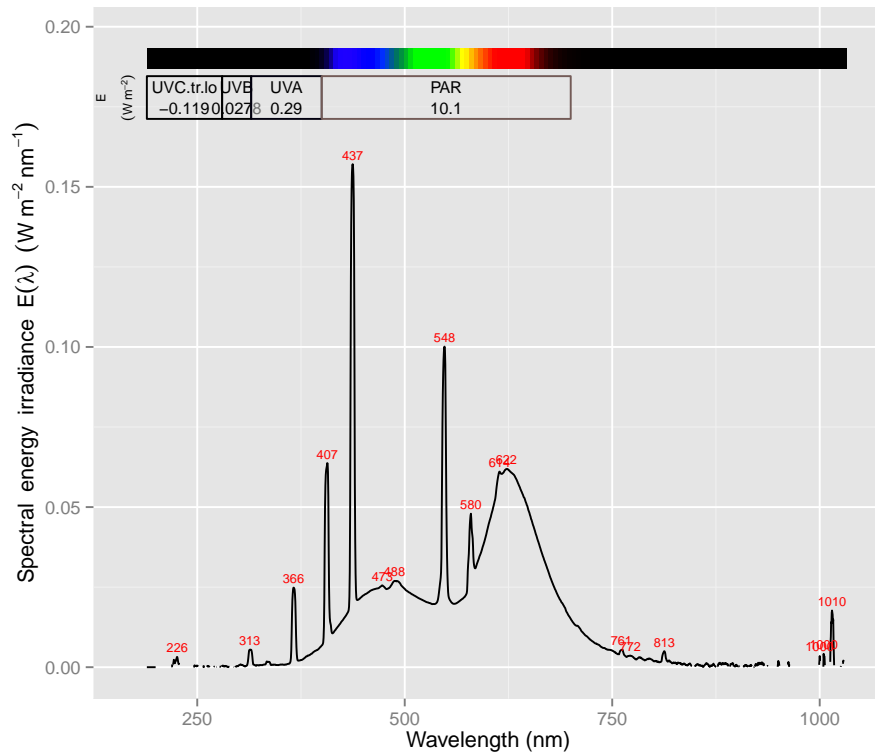

All in one statement.

```
plot(read_oोजaz_file(file = "spectrum.JazIrrad"))


## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning: Removed 1 rows containing missing values (geom_path).
```

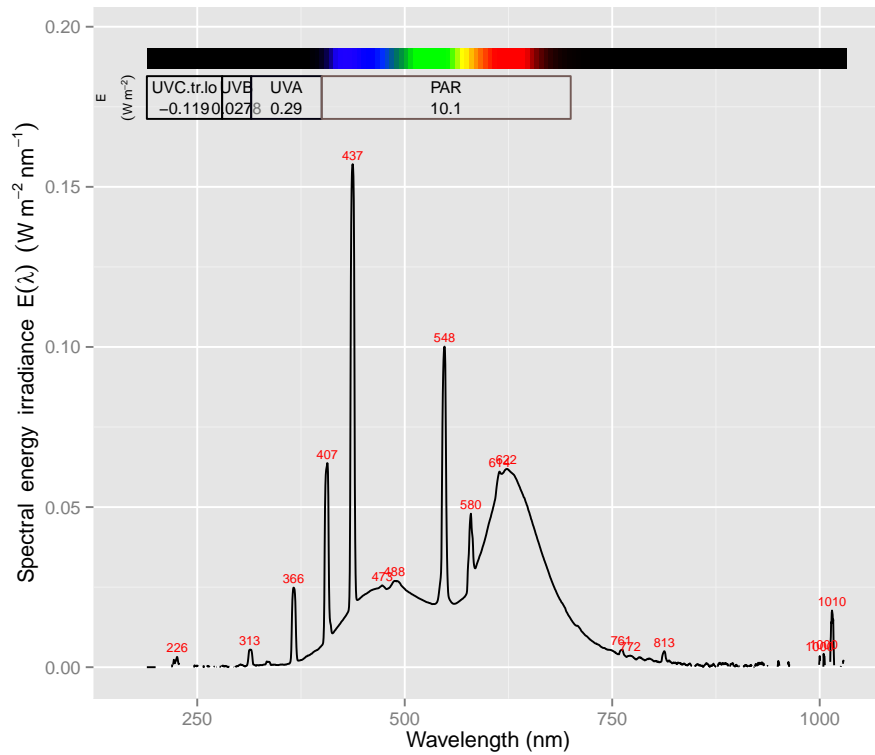

As above but limiting the wavelength range plotted.

```
plot(read_oobjaz_file(file = "spectrum.JazIrrad"), range = c(250,850))


## Warning in range.check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range.check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range.check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range.check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range.check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.00044
```

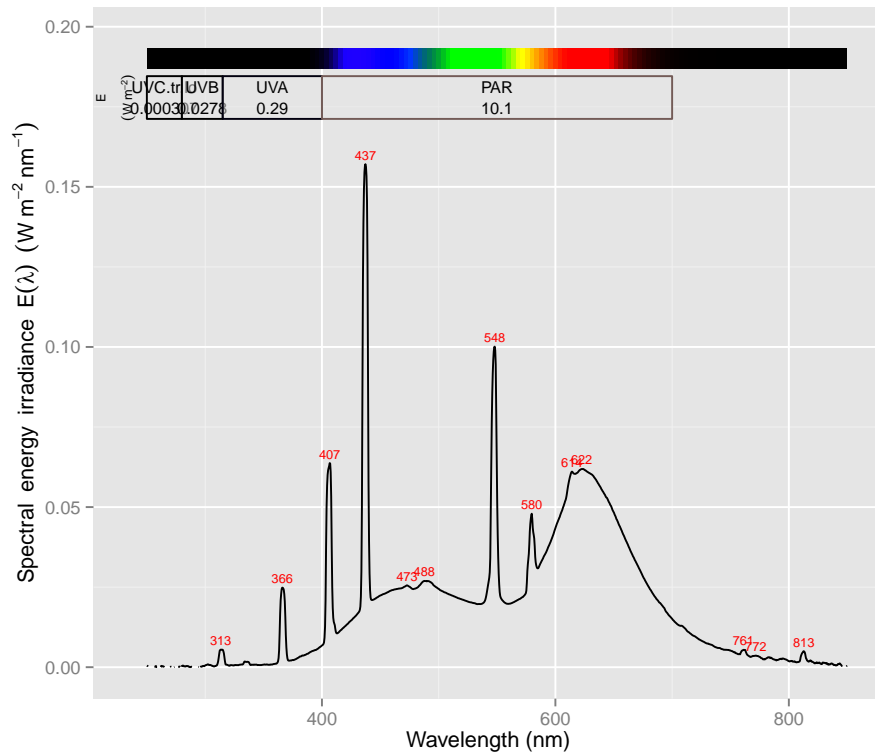

Adding our custom “adaptive” smoothing.


```
plot(smooth_spect(read_ojaz_file(file = "spectrum.JazIrrad")), range = c(315,900))

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in smooth_spect.source_spect(read_ojaz_file(file = "spectrum.JazIrrad")):
89 'bad' estimates in spectral irradiance
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.027
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.027
```

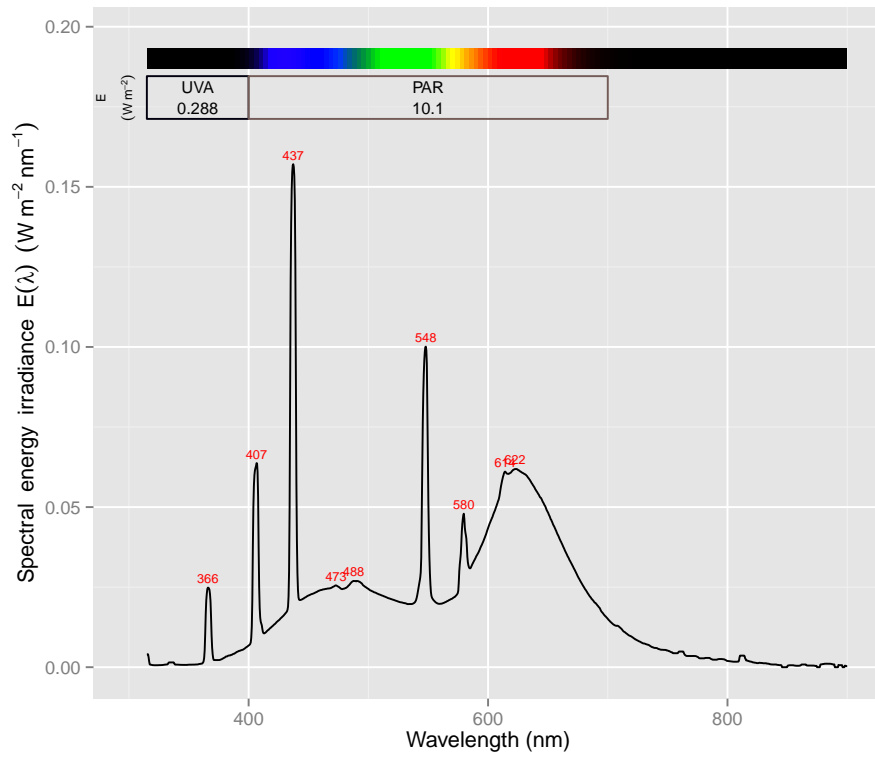


Now a file from an Ocean Optics' Q6500? spectrometer.

```
plot(read_ooss_file(file = "spectrum.SSIrrad"))
```

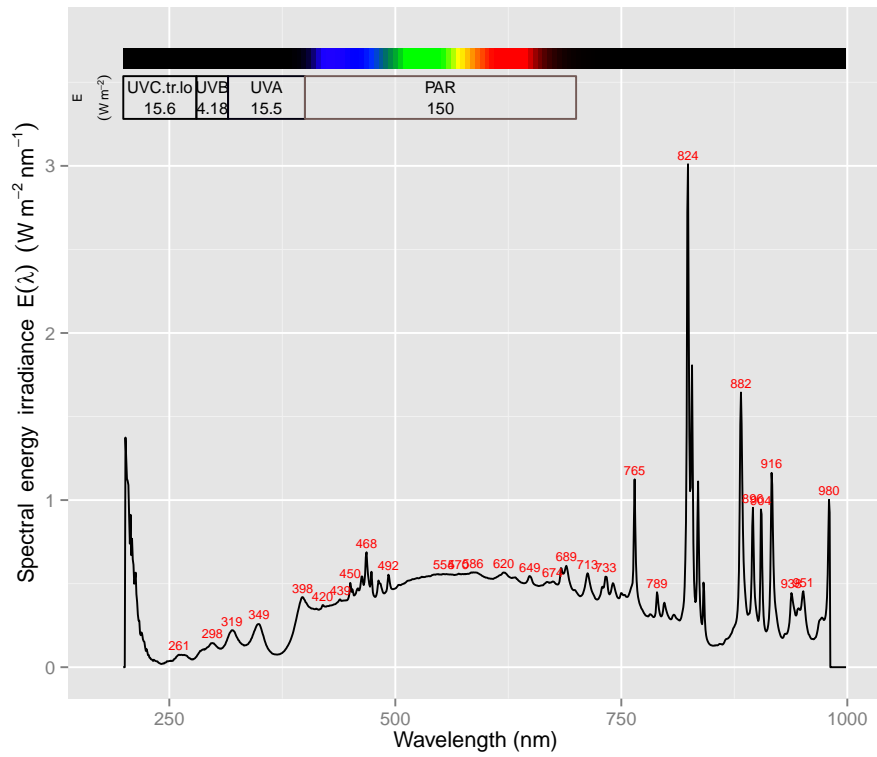


Macam's single column DTA files can also be imported.

```
plot(read_macam_file(file = "spectrum.DTA"))
```

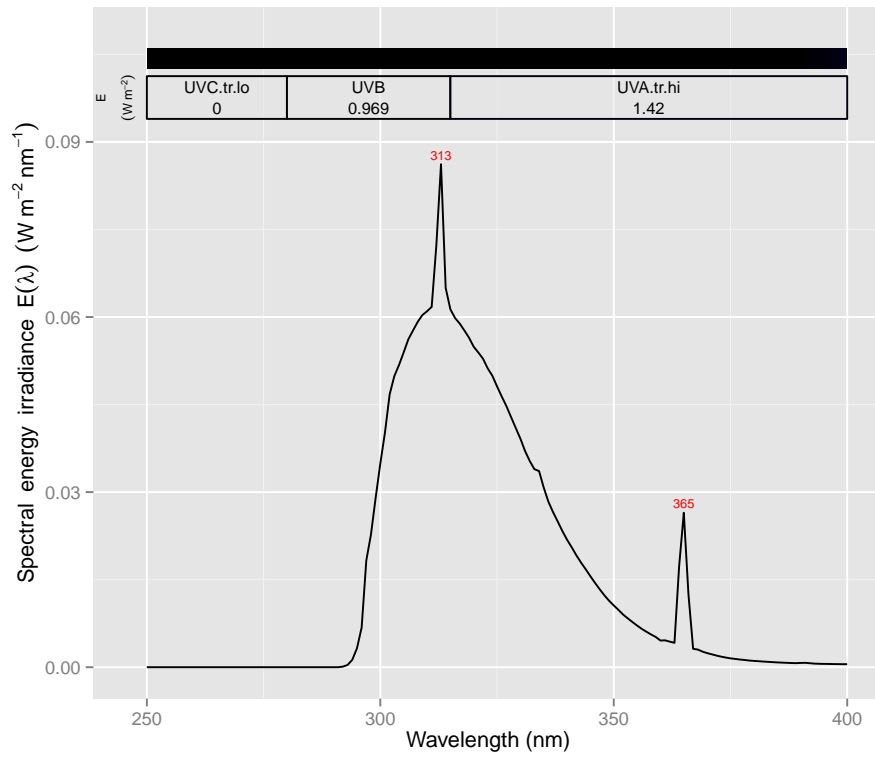

And a file generated by LI-COR's PC1800 program for the LI-1800 septro-radiometer.

```
plot(read_licor_file(file = "spectrum.PRN"))
```

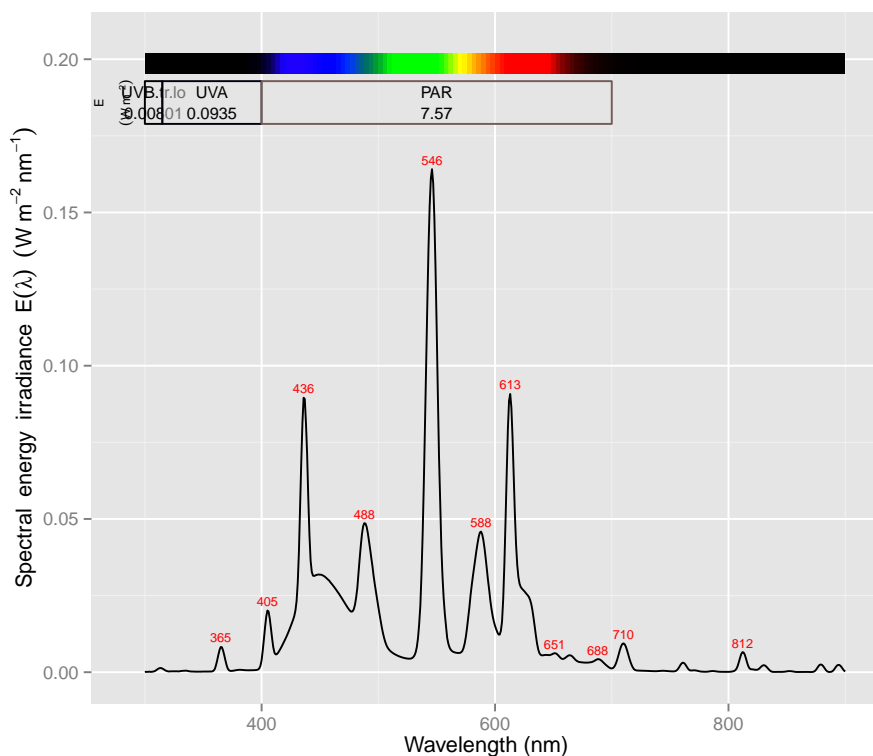In all cases as much information as possible is decoed, and the data file headers are preserved as comments in the source.spct objects.

```
cat(comment(jaz.spct))

## Ocean Optics:
## Jaz Absolute Irradiance File
## ++++++
## Date: Tue Feb 03 09:44:41 2015
## User: jaz
## Dark Spectrum Present: Yes
## Processed Spectrum Present: Yes
## Spectrometers: JAZA1065
## Integration Time (usec): 193000 (JAZA1065)
## Spectra Averaged: 3 (JAZA1065)
## Boxcar Smoothing: 5 (JAZA1065)
## Correct for Electrical Dark: Yes (JAZA1065)
## Strobe/Lamp Enabled: No (JAZA1065)
## Correct for Detector Non-linearity: Yes (JAZA1065)
## Correct for Stray Light: No (JAZA1065)
## Number of Pixels in Processed Spectrum: 2048
## Fiber (micron): 3900
## Collection Area: 0.119459
## Int. Sphere: No
```

A variable with the user supplied date and time data, or the date read from

the header (the text itself) not the file date as the file date may not reflect the creation date and time.

```
jaz.spct <- read_oojaz_file(file = "spectrum.JazIrrad", date = NULL)

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032

jaz.spct

## Object: local source_spct [2,047 x 3]
## Wavelength (nm): range 189.28485 to 1033.1483, step 0.357056 to 0.459564
##
##      w.length s.e.irrad date
## 1  189.2849      0    NA
## 2  189.7444      0    NA
## 3  190.2040      0    NA
## 4  190.6635      0    NA
## 5  191.1229      0    NA
## 6  191.5824      0    NA
## 7  192.0417      0    NA
## 8  192.5011      0    NA
## 9  192.9604      0    NA
## 10 193.4197      0    NA
##
```

The output from the TUV model can be imported either by editing it before import, or by making a simple edit to the output routine of TUV. This function is known to work with TUV version 5.0 output. The output from TUV can contain a variable number of spectra in “parallel” columns, which are *melted* into a single column, with a factor with letter as levels, a numeric variable with the zenith angle and a POSIXct column with times. A date needs to be always supplied as the output file from TUV has only time of day information.

```
tuv.spct <- read_tuv_file(file = "usrout.txt", date = ymd("2014-03-21"))

## Warning in check_generic_spct(x, multiple.wl = multiple.wl): 'w.length' values
## are not unique: 8 copies of each.

summary(subset(tuv.spct, spectrum == "A"))

## wavelength ranges from 280.5 to 761.5 nm
## largest wavelength step size is 1 nm
## spectral irradiance ranges from 3.041e-15 to 1.947 W m-2 nm-1
## energy irradiance is 636.7 W m-2

tuv.spct

## Object: local_source_spct [3,856 x 5]
## Wavelength (nm): range 280.5 to 761.5, step -481 to 1
##
```

```
##      w.length spectrum s.e.irrad angle      date
## 1      280.5      A 3.041e-15 1.829 2014-03-21 12:00:00
## 2      281.5      A 1.164e-13 1.829 2014-03-21 12:00:00
## 3      282.5      A 1.824e-12 1.829 2014-03-21 12:00:00
## 4      283.5      A 9.874e-12 1.829 2014-03-21 12:00:00
## 5      284.5      A 1.258e-10 1.829 2014-03-21 12:00:00
## 6      285.5      A 5.423e-10 1.829 2014-03-21 12:00:00
## 7      286.5      A 7.367e-09 1.829 2014-03-21 12:00:00
## 8      287.5      A 4.932e-08 1.829 2014-03-21 12:00:00
## 9      288.5      A 2.926e-07 1.829 2014-03-21 12:00:00
## 10     289.5      A 1.813e-06 1.829 2014-03-21 12:00:00
## ..      ...      ...      ...      ...      ...
```

It is possible to extract individual spectra with `subset`, or as done here plot them in different panels.

```
plot(tuv.spct, annotations = c("colour_guide")) +
  facet_wrap(~date, ncol = 2)
```

