photobiologyInOut Version 0.3.2 User Guide

Pedro J. Aphalo

October 2, 2015

Contents

1	Introduction	1
2	Evamples	-

1 Introduction

options(dplyr.print_max = 5)
options(dplyr.print_min = 3)

```
library(photobiologyWavebands)
library(photobiologygg)
library(photobiologyInOut)
library(lubridate)
library(ggplot2)
```

This package defines functions for importing spectral data from different instruments (Table 1) and simulation models (Table 2).

R function	Instrument	Program
read_ooss_file() read_oojaz_file() read_macam_file() read_licor_file()	Ocean Optics spectrometers Ocean Optics Jaz Macam LI-COR LI-1800	SpectraSuite instrument instrument PC1800 (MS-DOS)

Table 1: Functions for importing measured spectral emission data

R function	Simulation model	version
read_tuv_file() read_tuv_file()	TUV (Sasha Madronich) TUV (Sasha Madronich)	

Table 2: Functions for importing simulated spectral data from models

2 Examples

Reading a file geenrated by Ocean Optics' Jaz spectrometer.

```
jaz.spct <- read_oojaz_file(file = "spectrum.JazIrrad")

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032</pre>
```

Photon and energy irradiances.

```
q_irrad(jaz.spct, PAR()) * 1e6 # umol m-2 s-1
##
        PAR
## 47.43464
## attr(,"time.unit")
## [1] "second"
## attr(, "radiation.unit")
## [1] "photon irradiance total"
e_irrad(jaz.spct, PAR())
                              # W m-2
##
        PAR
## 10.12685
## attr(,"time.unit")
## [1] "second"
## attr(,"radiation.unit")
## [1] "energy irradiance total"
```


Plotting the spectrum.

```
plot(jaz.spct)

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032


## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032

## Warning: Removed 1 rows containing missing values (geom_path).
```

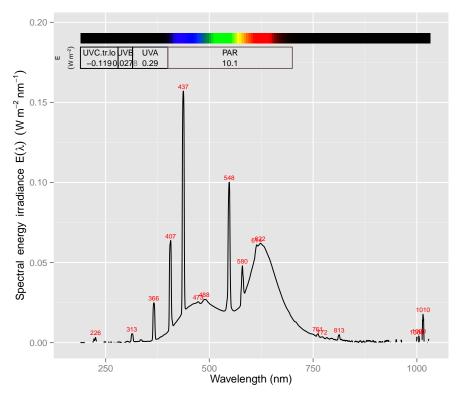


All in one statement.

```
plot(read_oojaz_file(file = "spectrum.JazIrrad"))
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning: Removed 1 rows containing missing values (geom_path).
```

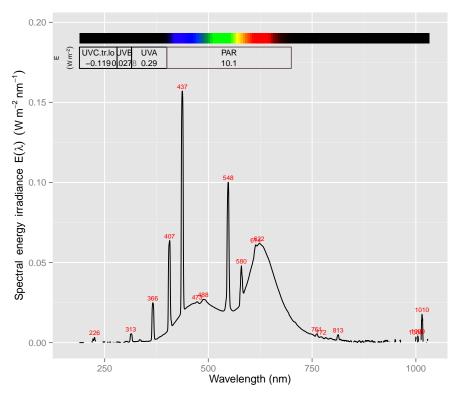

As above but limiting the wavelength range plotted.

```
plot(read_oojaz_file(file = "spectrum.JazIrrad"), range = c(250,850))
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimum s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range):
                                                            Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range):
                                                            Negative spectral energy
irradiance values; minimum s.e.irrad = -0.00044
## Warning in range_check(x, strict.range = strict.range):
                                                            Negative spectral energy
irradiance values; minimun s.e.irrad = -0.00044
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.00044
```

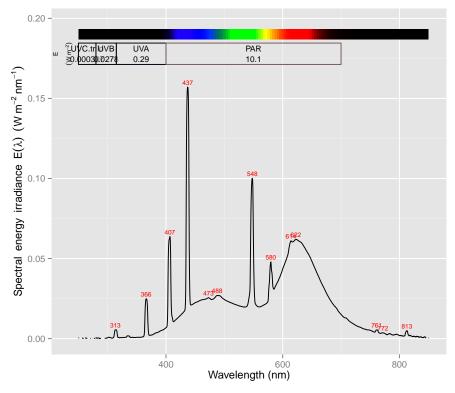

Adding our custom "adaptive" smoothing.

```
plot(smooth_spct(read_oojaz_file(file = "spectrum.JazIrrad")), range = c(315,900))


## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032

## Warning in smooth_spct.source_spct(read_oojaz_file(file = "spectrum.JazIrrad")):
89 'bad' estimates in spectral irradiance


## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.027

## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.027
```

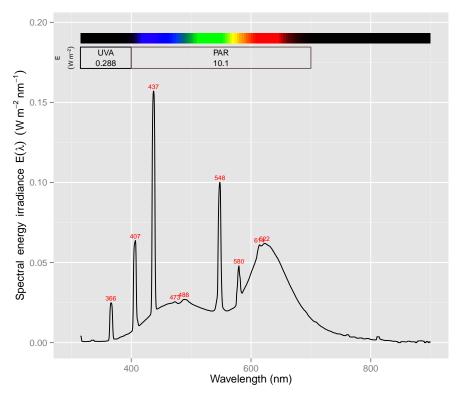


Now a file from an Ocean Optics' Q6500? spectrometer.

```
plot(read_ooss_file(file = "spectrum.SSIrrad"))
```

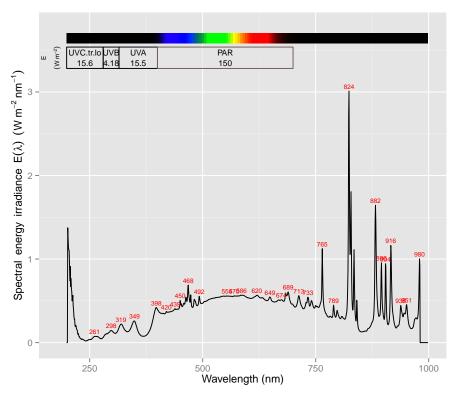


Macam's single column DTA files can also be imported.

```
plot(read_macam_file(file = "spectrum.DTA"))
```

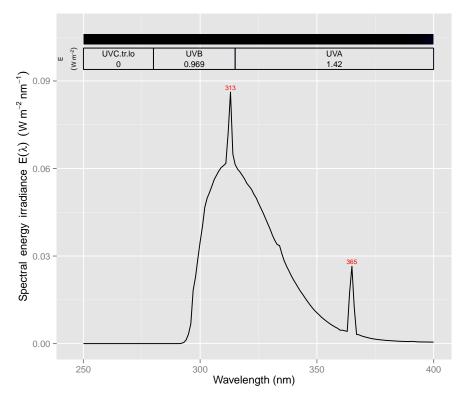

And a file generated by LI-COR's PC1800 program for the LI-1800 sepctroradiometer.

```
plot(read_licor_file(file = "spectrum.PRN"))
```

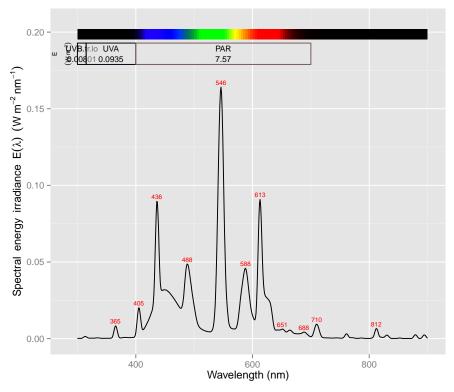

In all cases as much information as possible is decoed, and the data file headers are preserved as comments in the source.spct objects.

```
cat(comment(jaz.spct))
## Ocean Optics:
## Jaz Absolute Irradiance File
## Date: Tue Feb 03 09:44:41 2015
## User: jaz
## Dark Spectrum Present: Yes
## Processed Spectrum Present: Yes
## Spectrometers: JAZA1065
## Integration Time (usec): 193000 (JAZA1065)
## Spectra Averaged: 3 (JAZA1065)
## Boxcar Smoothing: 5 (JAZA1065)
## Correct for Electrical Dark: Yes (JAZA1065)
## Strobe/Lamp Enabled: No (JAZA1065)
## Correct for Detector Non-linearity: Yes (JAZA1065)
## Correct for Stray Light: No (JAZA1065)
## Number of Pixels in Processed Spectrum: 2048
## Fiber (micron): 3900
## Collection Area: 0.119459
## Int. Sphere: No
```

A variable with the user supplied date and time data, or the date read from

the header (the text itself) not the file date as the file date may now reflect the creation date and time.


```
jaz.spct <- read_oojaz_file(file = "spectrum.JazIrrad", date = NULL)</pre>
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Warning in range_check(x, strict.range = strict.range): Negative spectral energy
irradiance values; minimun s.e.irrad = -0.032
## Object: source_spct [2,047 x 3]
## Wavelength (nm): range 189.28485 to 1033.1483, step 0.357056 to 0.459564
## Time unit: 1s
##
##
      w.length s.e.irrad
                                        date
##
       (dbl) (dbl)
                                      (time)
                  0 2015-02-03 09:44:41
0 2015-02-03 09:44:41
## 1 189.2849
## 2 189.7444
## 3 190.2040
                    0 2015-02-03 09:44:41
```

The output from the TUV model can be imported either by editing it before import, or by making a simple edit to the output routine of TUV. This function is known to work with TUV version 5.0 output. The output from TUV can contain a variable number of spectra in "parallel" columns, which are *melted* into a single column, with a factor with letter as levels, a numeric variable with the zenith angle and a POSIXct column with times. A date needs to be always supplied as the output file from TUV has only time of day information.

```
tuv.spct <- read_tuv_file(file = "usrout.txt", date = ymd("2014-03-21"))</pre>
summary(subset(tuv.spct, spct.idx == "A"))
## wavelength ranges from 280.5 to 761.5 nm
## largest wavelength step size is 1 nm
## spectral irradiance ranges from 3.041e-15 to 1.947 W m-2 nm-1
## energy irradiance is 636.7 \text{ W m-}2
tuv.spct
## Object: source_spct [3,856 x 5]
## Wavelength (nm): range 280.5 to 761.5, step -481 to 1
## Time unit: 1s
##
##
      w.length spct.idx s.e.irrad angle
        (dbl) (fctr) (dbl) (dbl)
##
                                                    (time)
                 A 3.041e-15 1.829 2014-03-21 12:00:00
## 1
         280.5
## 2
                    A 1.164e-13 1.829 2014-03-21 12:00:00
        281.5
## 3
        282.5
                    A 1.824e-12 1.829 2014-03-21 12:00:00
```

It is possible to extract individual spectra with subset, or as done here plot them in different panels.

```
plot(tuv.spct, annotations = c("colour.guide")) +
  facet_wrap(~date, ncol = 2)
```

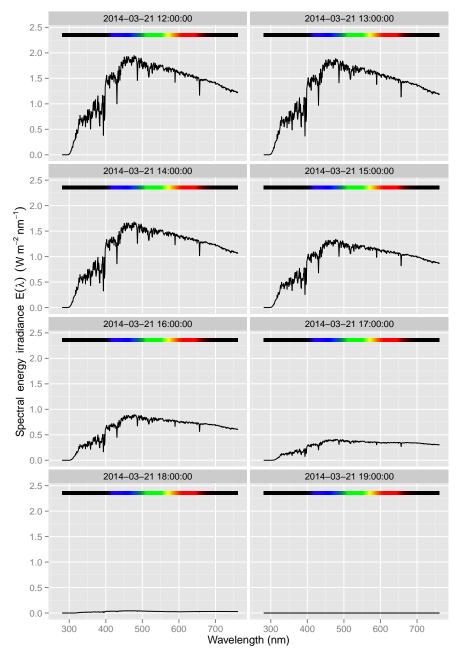

The output is a single source_spct object that can be easily converted into a source_mspct object containing the individual spectra as members of the

collection.

```
tuv.mspct <- subset2mspct(tuv.spct)</pre>
tuv.mspct
## Object: source_mspct [8 x 1]
## --- Member: A ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
     w.length s.e.irrad angle
                                            date
        (dbl) (dbl) (dbl)
##
## 1
        280.5 3.041e-15 1.829 2014-03-21 12:00:00
## 2
      281.5 1.164e-13 1.829 2014-03-21 12:00:00
      282.5 1.824e-12 1.829 2014-03-21 12:00:00
## 3
## .. ... ... ... ... ... ## --- Member: B ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
     w.length s.e.irrad angle
                                             date
        (dbl) (dbl) (dbl)
## 1
        280.5 1.314e-15 13.198 2014-03-21 13:00:00
## 2
        281.5 5.415e-14 13.198 2014-03-21 13:00:00
## 3
        282.5 9.039e-13 13.198 2014-03-21 13:00:00
## ..
         ... ...
## --- Member: C ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
    w.length s.e.irrad angle
##
       (dbl) (dbl) (dbl)
                                          (time)
        280.5 4.521e-17 28.2 2014-03-21 14:00:00 281.5 2.510e-15 28.2 2014-03-21 14:00:00
## 1
## 2
      282.5 5.413e-14 28.2 2014-03-21 14:00:00
## 3
## ..
        ...
## --- Member: D ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1 \,
## Time unit: 1s
##
##
     w.length s.e.irrad angle
                                             date
##
        (dbl) (dbl) (dbl)
        280.5 3.075e-20 43.202 2014-03-21 15:00:00
## 1
        281.5 3.273e-18 43.202 2014-03-21 15:00:00
## 2
## 3 282.5 1.234e-16 43.202 2014-03-21 15:00:00
## ..
       ... ...
## --- Member: E ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
     w.length s.e.irrad angle
## (dbl) (dbl) (dbl)
                                            (time)
```

```
## 1 280.5 2.253e-26 58.205 2014-03-21 16:00:00
         281.5 7.751e-24 58.205 2014-03-21 16:00:00
## 2
## 3
         282.5 8.148e-22 58.205 2014-03-21 16:00:00
## .. ... ...
## --- Member: F ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1 \,
## Time unit: 1s
##
##
    w.length s.e.irrad angle
##
       (dbl) (dbl) (dbl)
## 1
         280.5 1.929e-27 73.208 2014-03-21 17:00:00
## 2
         281.5 5.710e-25 73.208 2014-03-21 17:00:00
## 3
       282.5 5.202e-23 73.208 2014-03-21 17:00:00
        ...
## --- Member: G ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
      w.length s.e.irrad angle
                                                date
##
        (dbl) (dbl) (dbl)
## 1
        280.5 4.721e-28 88.211 2014-03-21 18:00:00
## 2
        281.5 1.385e-25 88.211 2014-03-21 18:00:00
## 3 282.5 1.250e-23 88.211 2014-03-21 18:00:00
## ..
## --- Member: H ---
## Object: source_spct [482 x 4]
## Wavelength (nm): range 280.5 to 761.5, step 1
## Time unit: 1s
##
##
     w.length s.e.irrad angle
  (dbl) (dbl) (dbl)
                                                 date
##
      280.5 0 103.213 2014-03-21 19:00:00

281.5 0 103.213 2014-03-21 19:00:00

282.5 0 103.213 2014-03-21 19:00:00

...
                                               (time)
## 1
## 2
## 3
## ..
                    ...
## --- END ---
```

With the default of lubridate::today() for date times are 'mapped' to the current local date using the time zone of the computer as visible to R.

```
tuv_nd.spct <- read_tuv_file(file = "usrout.txt")</pre>
tuv_nd.spct
## Object: source_spct [3,856 x 5]
## Wavelength (nm): range 280.5 to 761.5, step -481 to 1
## Time unit: 1s
##
##
      w.length spct.idx s.e.irrad angle
       (dbl) (fctr) (dbl) (dbl)
                                                         (time)
         280.5 A 3.041e-15 1.829 2015-10-02 12:00:00
## 1
         281.5 A 1.164e-13 1.829
282.5 A 1.824e-12 1.829
... ...
                      A 1.164e-13 1.829 2015-10-02 12:00:00
A 1.824e-12 1.829 2015-10-02 12:00:00
## 2
## 3
## ..
```